

SYST'AM® P914L / UNIVERSAL CUSHION


POSTURAL AID AND PRESSURE ULCER PREVENTION

MATERIALS

80
kg/m³



Visco foam

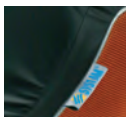

Removable
POLYMAILLE®
integral cover



Foam maintenance:



(Do not immerse in water)

Cleaning POLYMAILLE® cover:



INDICATIONS


- Universal positioning help.



VISCOELASTIC FOAM WITH SHAPE MEMORY


- Defined body moulding and increased body surface area in contact with the foam.
- More effective distribution of pressure: reduction of transcutaneous pressure on highrisk zones, facilitated blood flow.
- Increased comfort and more effective maintenance of patient position.

AVAILABLE VERSIONS

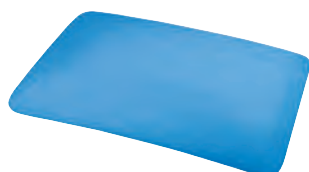

P914L / UNIVERSAL CUSHION

- This is the perfect supporting aid: it helps support the patient, relieves the pressure on the heels, protects the internal lateral knee areas, ankles...
- This cushion can also be used as postural help in Rest chair, wheelchair and shell seat.
- Facilitates positioning for care staff.
- Used by itself or with other devices.

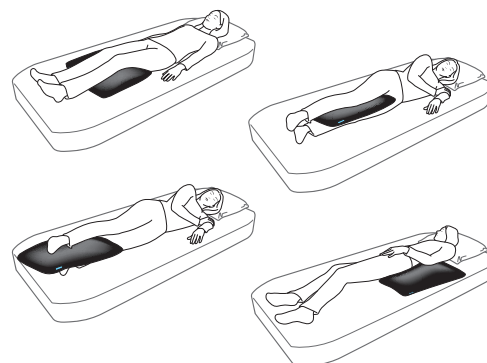

FEATURES OF THE COVERS

FABRIC COATED WITH BI-STRETCH POLYURETHANE

- Reduces friction and shear effects.
- Supple and soft to the touch (comfortable).
- Favours the exchange of gases (steam, sweat):
 - fights against maceration.
- Impermeable material:
 - better hygiene,
 - longer support system lifespan.
- Washable at 90°C, can be decontaminated using cold sprays.
- Treated to resist fire.
- In multi-patient use, it is preferable to buy one new cover per patient.



DESIGNATION	ITEM CODE	SIZE L x W x H (cm / inches)
P914L / Universal cushion	P914L1HW	60 x 40 x 13 cm / 24 x 15,7 x 5"

